



Canis Major[®] Exam Lift Table

For models:
CMWO (without scale)
CMWS (with scale)



User's Guide

Product Information

(The information below is required when calling for service)


Dealer [name / phone]:

Date of Purchase:

Model / Serial Number:

Midmark Authorized Service Company [name / phone]:

Model / Serial Number Location

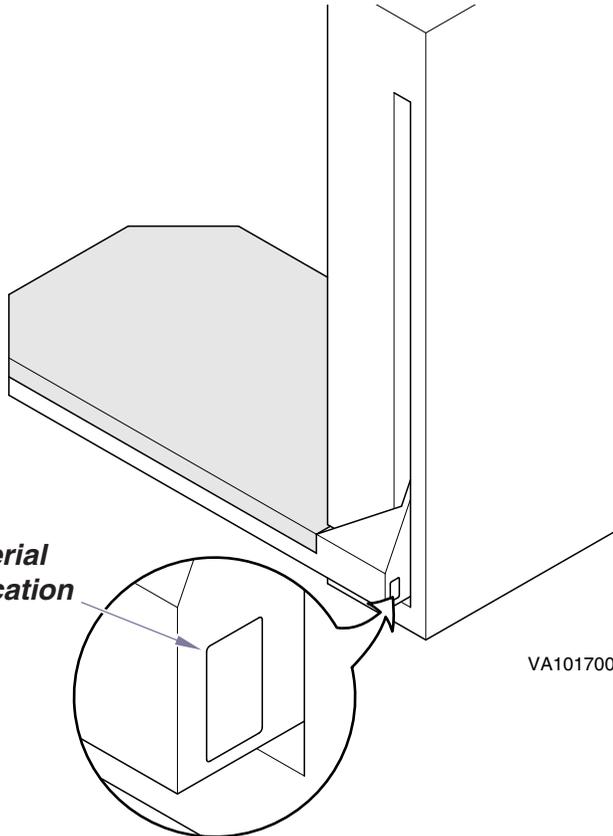

VA101700i

Table of Contents

Important Information

| | |
|---|---|
| Safety Symbols | 4 |
| Transportation / Storage Conditions | 4 |
| Intended Use | 4 |
| Electrical Precautions | 4 |

Operation

| | |
|---------------------------------------|---|
| Table Up / Down (Hand Switches) | 5 |
| Table Up / Down (Foot Control) | 5 |
| Scale (Optional) | 6 |

Maintenance

| | |
|--------------------------------|---|
| Cleaning | 8 |
| Preventative Maintenance | 8 |
| Troubleshooting | 9 |
| Calling For Service | 9 |

Specifications

| | |
|--|----|
| Weights, Dimensions, Electrical Specifications | 10 |
| Range of Motion | 10 |

Warranty Information

| | |
|------------------------|----|
| Limited Warranty | 11 |
|------------------------|----|

Important Information

Safety Symbols



DANGER

Indicates an imminently hazardous situation which will result in serious or fatal injury if not avoided. This symbol is used only in the most extreme conditions.



WARNING

Indicates a hazardous situation which could result in serious injury if not avoided.



CAUTION

Indicates a potentially hazardous situation which could result in minor or moderate injury if not avoided. It may also be used to alert against unsafe practices.



Equipment Alert

Indicates a potentially hazardous situation which could result in equipment damage if not avoided.

Note

Amplifies a procedure, practice, or condition.



Dangerous voltage / shock hazard



Consult User Guide for important information.

Transportation / Storage Conditions

Ambient Temperature Range: 5°C to +38°C (+41°F to 100°F)

Relative Humidity 10% to 90% (non-condensing)

Atmospheric Pressure 500hPa to 1060hPa (0.49atm to 1.05atm)

Intended Use

This product is intended to be used as a table to provide positioning, support and cleaning of animals during examination or minor surgery procedures conducted by qualified professionals.

Electrical Precautions



WARNING

When using high frequency surgical devices or endocardial catheters:

- Use nonconductive material to insulate patient from metal portions of table.
- Consult operating instructions for the device before use.

Failure to comply may result in electric shock or burns to patient.

Operation



CAUTION


If the down function stops during operation, immediately release Down control and remove patient from table. Stand clear of tabletop and use the Up control to raise the table and remove obstruction. Ensure area around table is clear and use the Down control to lower table to the desired height. Failure to comply with these instructions may result in injury to the patient, staff and/or equipment.

Table Up / Down (Hand Switches)

To move the table top UP or DOWN...
Press and hold either switch. Press the upper portion to raise table or lower portion to lower.

Note

Movement automatically stops at the max. / min. height. This table is also equipped with a safety feature that disables the down function when something is encountered under the table.

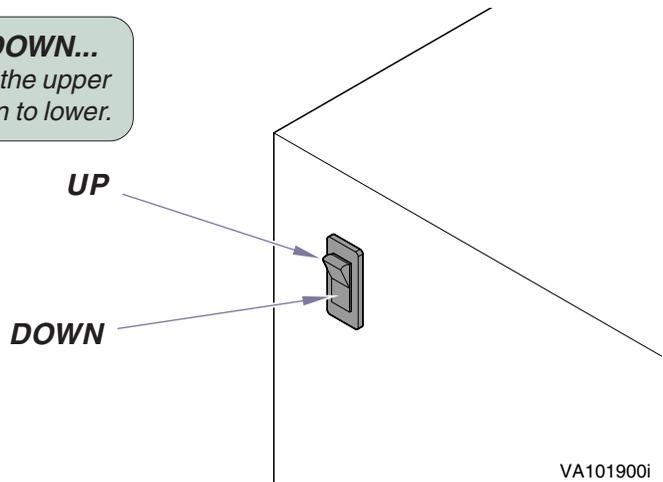


VA101900i

Table Up / Down (Foot Control)

To move the table top UP or DOWN...
Press and hold the appropriate switch.

Note

Movement automatically stops at the max. / min. height. This table is also equipped with a safety feature that disables the down function when something is encountered under the table.

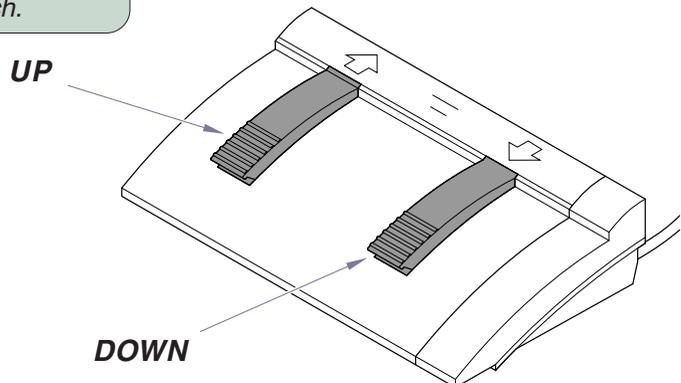


VA102000i

Operation

Scale (optional)

The 300 pound capacity platform scale utilizes four (4) built in load sensors. These electronic sensors are accurate within 5 grams, and are protected against impacts, moisture and liquid spills. The scale has been calibrated and all parameters are preset so it should never need calibration.

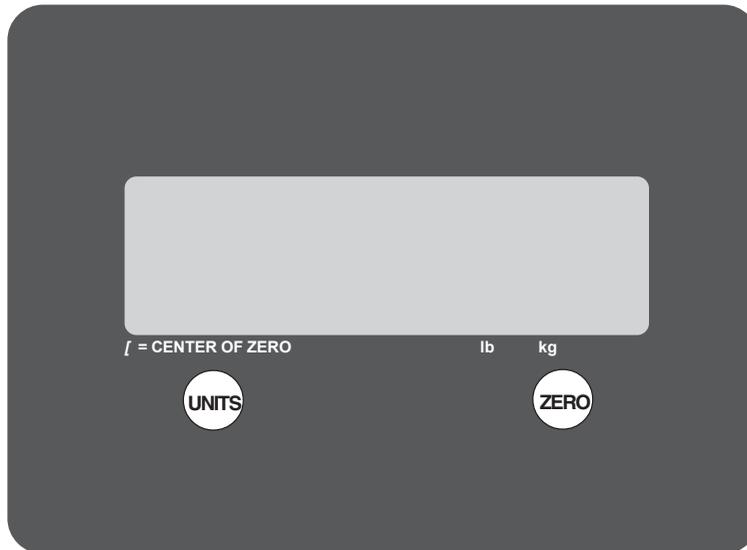

Equipment Alerts

- Prevent flammables and liquids from entering the scale head.
- Allow adequate clearance on all sides of scale platform to ensure accuracy.
- Do not drop large loads onto the scale platform.
- Never exceed the rated capacity of the scale.
- Do not operate in close proximity to RF transmitters.

There are two (2) controls located on the display:

UNITS Toggles the displayed weight between pounds (lbs) and kilograms (kg).
ZERO Resets the display to indicate zero (0.0).

Note

The scale will not respond to the ZERO control under these conditions:

Motion on Scale

An under load condition (less than 0.1lbs or 50g)

An over load condition (more than 300lbs or 1kg)

Operation

Scale (continued)

Weighing an animal without a cage:

- Zero scale (0.0), place animal on center of platform and read the weight displayed.

Note

A blank scale while weighing means excess motion is detected and weight averaging is in process. Once completed, the scale will display and HOLD the recorded weight. Adding additional weight while in the HOLD condition will NOT change the display.

Weighing an animal with a cage:

- Place empty cage on center of platform, then zero scale (0.0).
- Remove cage from platform; scale will read negative weight (-_._).
- Place cage with animal on center of platform and read/record weight.
- To continue to weigh animals in same sized cage, place next cage with animal on center of platform and read/record weigh; do not zero scale.

Troubleshooting:

| <u>Symptom or Error</u> | <u>Possible Problem (and/or) Remedy</u> |
|---------------------------------------|---|
| No Display | Power is OFF. / Table is unplugged from the power source. |
| Incorrect weight | Interference with scale platform. / Ensure there is nothing in contact with the sides of the scale platform. Ensure the correct UNIT (of measure) is selected. |
| Display segment (s) missing/illegible | Low or intermittent voltage. / Unplug for one (1) minute to reset system then reconnect to original power source. Disconnect power cord and use different power source. |
| ZERO key does not reset to [0.0] | Motion or activity on the scale. / Remove source of motion. Ensure platform is level and supported evenly at corners. Weight on scale may exceed the zero limit. Unplug for one (1) minute to reset system. |
| Display “locked & will not change | System Error. / Unplug the unit for one (1) minute to reset system. then reconnect to original power source. |
| Displays “HiLoAd” | Weight on scale exceeds capacity. / Remove excess load. |
| Displays “HiLdX” | The load placed on the scale is too heavy on a corner [shown as X]. / Move load to center of the scale platform. |
| Displays “LoLdX” | The scale is under loaded / Ensure that tabletop cover pins are properly seated on load cells. |

If the problem cannot be resolved using the above corrective actions, please contact Midmark at 1-800-Midmark (1-800-643-6275) or 937-526-3662.

Maintenance

Cleaning

Tabletop Mat:

The tabletop mat should be cleaned between patients with a virucidal agent or mild disinfectant spray and towel or cloth dried. The mat can be moved to a tub or sink for a more thorough cleaning.



Equipment Alert

Silver nitrate and betadine can stain the mat. Certain pharmaceutical products such as “ADVANTAGE” will react and cause bubbles on the mat surface. Caution must be used when selecting a cleaning agent.

Painted Surfaces:

Clean the table weekly, wiping the painted metal and plastic surfaces with a clean, soft cloth, and mild cleaners. Periodic applications of common furniture wax will ease cleaning and maintain the finished luster of the table’s painted surfaces.

Stainless Steel Surfaces:

For best results, stainless steel should be cleaned weekly or whenever films or deposits become apparent. Use a mild cleaning agent and whenever possible rinse thoroughly and towel dry. More stubborn deposits or tightly adhering debris may require vigorous rubbing, or brushing. If a commercial cleaner needs to be used, make sure it is approved for use on stainless steel surfaces.



Equipment Alert

Always use nylon or plastic bristle scrubbers. The use of steel wool or steel brushes may cause steel particles to become embedded in the surfaces and will cause rust.

Avoid prolonged surface exposure to chlorides, bromides, iodides and thiocyanates. Salty solutions allowed to evaporate on stainless steel surfaces may cause pitting, corrosion or metal discoloration.

Preventive Maintenance

Periodically inspect the following areas:

- Power cord should be free of cuts or other visible damage.
- All fasteners should be present and tightened securely.
- All mechanical functions should operate properly.

Have an authorized service technician inspect your table on a regular basis.

Maintenance

Troubleshooting

If the table will not go up, check the following:

- Unit is plugged into a wall outlet
- Power to the outlet
- No visible damage to wiring
- The table is at its maximum height
- Foot control is fully inserted and making proper connections (try either hand toggle switches)
- Toggle switch is bad (try using switch on other side or foot control)

If the table will not go down, check the following:

- Unit is plugged into a wall outlet
- Power to the outlet
- No visible damage to wiring
- The table is at its minimum height
- The safety pan switch is activated by an obstacle beneath the table
- The safety pan switch is out of adjustment
- Foot control is fully inserted and making proper connections (try either hand toggle switches)
- Toggle switch is bad (try using switch on other side or foot control)

If the table won't stop when encountering an obstacle beneath the table, check:

- The safety pan bound or switch is out of adjustment

If the scale does not operate, check the following:

- Refer to scale section elsewhere in this manual

Calling For Service

Note

Model / Serial Number information is required when calling for service.

If service is required, contact your Midmark dealer.

To contact Midmark directly:

1-800-Midmark (1-800-643-6275) or 937-526-3662
8:00 am until 5:00 pm. Monday through Friday (EST)
[excluding standard U.S. holidays]

Specifications

Weights, Dimensions, Electrical Specifications

Maximum Patient Weight: 300 lbs (136 kg)

Weight of Table:

45" tabletop (w/o scale) [unpacked] 315 lbs (142 kg)

45" tabletop (w/ scale) [unpacked] 320 lbs (145 kg)

45" tabletop (w/o scale) [with packaging & skid] 372 lbs (169 kg)

45" tabletop (w/ scale) [with packaging & skid] 377 lbs (171 kg)

Range of Motion & Dimensions: (see illustration below)

Electrical Requirements: Dedicated Circuit, 115 VAC, 60Hz, 10A



Power Cord: 4 ft. (121 cm) long

Duty Cycle

(actuator run time): Intermittent Operation

[1 minute ON / 9 minutes OFF]

Range of Motion & Dimensions

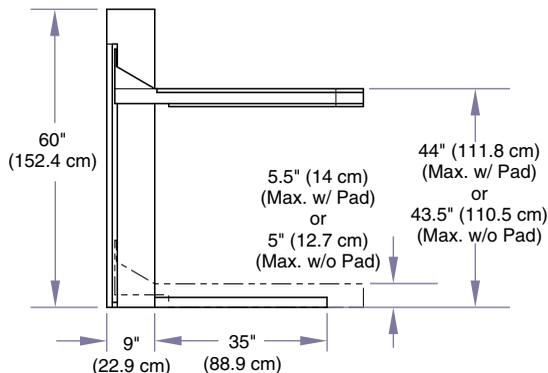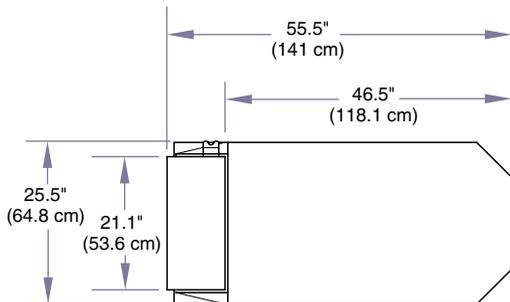

VA1700i

Warranty Information

SCOPE OF WARRANTY Midmark Corporation (“Midmark”) warrants to the original retail purchaser that it will repair or replace components of the animal health products manufactured by Midmark (except for products and components not warranted under “Exclusions”) that are defective in material or workmanship under normal use and service. The sole remedy under this limited warranty is the repair or replacement, at Midmark’s option, of the applicable products or components. This limited warranty shall only apply to defects that: (i) are reported to Midmark within the applicable warranty period; and (ii) are determined to exist upon examination by Midmark. This limited warranty extends only to the original retail purchaser of a product, and is not transferable or assignable.

APPLICABLE WARRANTY PERIOD The applicable warranty period, measured from the date of delivery to the original user, shall be one (1) year for all warranted products and components.

OBTAINING WARRANTY SERVICE Warranty service must be obtained through either Midmark or an authorized dealer in the Midmark product line for which warranty service is requested. Midmark may be contacted for warranty service inquiries or issues via email at www.midmark.com, by mail to Midmark Corporation, 60 Vista Drive, Versailles, Ohio 45380, or by phone at: 1-800-MIDMARK.

It is the retail purchaser’s obligation to arrange for delivery of a product to Midmark or one of its authorized dealers for warranty service, which delivery shall be at retail purchaser’s expense. It is also the retail purchaser’s obligation to comply with the warranty service instructions provided either by Midmark or its authorized dealer. The retail purchaser must provide Midmark with completed warranty registration information within thirty (30) days after purchase in order to obtain the benefits of this limited warranty.

EXCLUSIONS This limited warranty does not cover and Midmark shall not be liable for the following:

- (1) defects, damage or other conditions caused, in whole or in part, by misuse, abuse, negligence, alteration, accident (including animal acts of any kind), freight damage, tampering or failure to seek and obtain repair or replacement in a timely manner;
- (2) matching of color, grain or texture except to commercially acceptable standards;
- (3) changes in color caused by natural or artificial light;
- (4) products which are not installed, used and properly cleaned and maintained as required in the installation and operation manuals for the applicable product (imaging product must be installed by a certified Midmark installer);
- (5) products considered to be of a consumable nature;
- (6) accessories or parts not manufactured by Midmark;
- (7) specially manufactured products;
- (8) charges by anyone (including Midmark’s authorized dealers) for adjustments, repairs, replacement parts, installation or other work performed upon or in connection with such products which are not expressly authorized in writing in advance by Midmark;
- (9) costs and expenses of routine maintenance and cleaning;
- (10) all sinks, faucets and plumbing accessories;
- (11) representations and warranties made by any person or entity other than Midmark; and
- (12) with respect to software that is a product or a component thereof, that the software will be error free, can be used without problems or interruptions, or will be free from vulnerability to intrusion or attack by viruses or other methods.

EXCLUSIVE REMEDY; CONSEQUENTIAL DAMAGES DISCLAIMER Midmark’s only obligation under this limited warranty is the repair or replacement of defective parts. Midmark shall not be liable for and hereby disclaims any direct, special, indirect, incidental, exemplary or consequential damages or delays including, but not limited to, damages for loss of profits or income, loss of use, downtime, cover, and employee or independent contractor wages, payments and benefits. This disclaimer shall survive any failure or asserted failure of the essential purpose of this limited warranty or its remedies specified herein.

NO AUTHORIZATION No person or firm is authorized to create or approve for Midmark any other obligation or liability in connection with Midmark products.

WARRANTY DISCLAIMER: THIS LIMITED WARRANTY IS MIDMARK’S ONLY WARRANTY AND IS IN LIEU OF ALL OTHER WARRANTIES, EXPRESS OR IMPLIED. MIDMARK MAKES NO IMPLIED WARRANTIES OF ANY KIND INCLUDING ANY IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE.

STATUTE OF LIMITATIONS No action may be brought against Midmark for breach of this limited warranty, an implied warranty, if any, or for any other claim arising out of or relating to the products, more than ninety (90) days following expiration of the warranty period. In the event multiple warranty periods exist with respect to a product, the ninety (90) day period provided for herein shall begin to run from expiration of the warranty period for the component to which the claim relates.

SEVERABILITY In the event any provision of this limited warranty is determined to be invalid or otherwise unenforceable: (i) the provision shall be enforced in a manner that closest holds to its intent, while at the same time curing the invalidity or unenforceability, and (ii) the balance of this limited warranty, being severable, shall not be affected in the event the provision that is invalid or unenforceable cannot be enforced in any respect.

Midmark Corporation
60 Vista Drive
Versailles, OH 45380-0286
1-800-643-6275
1-937-526-3662
www.midmark.com



Because we care.