
Ritter®

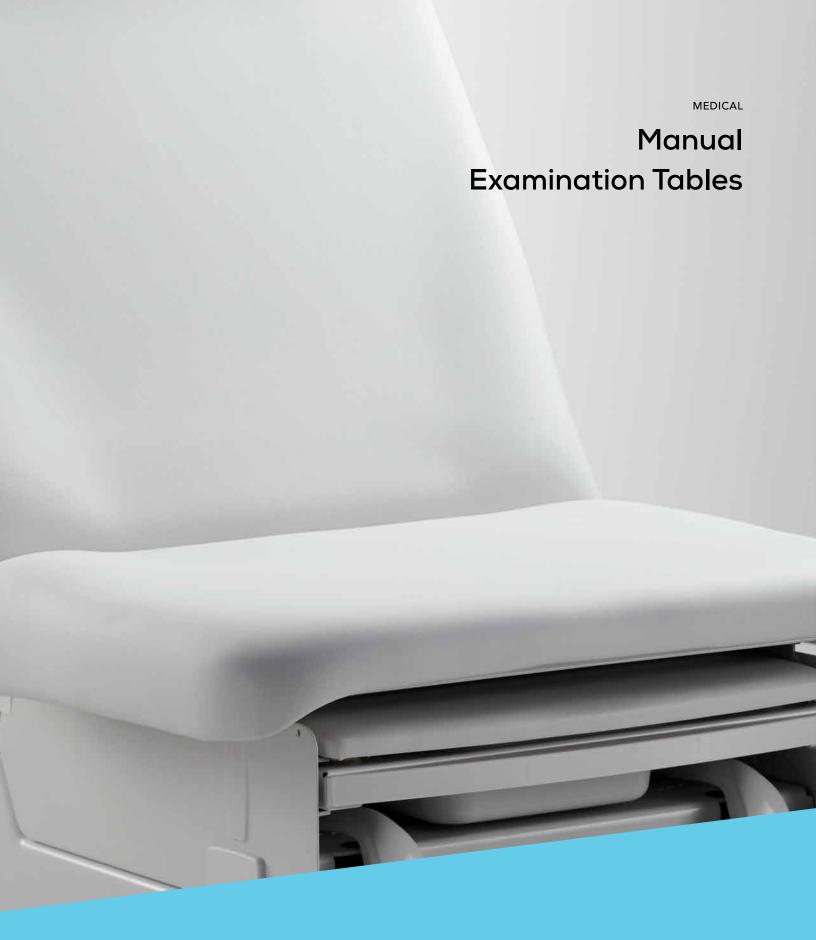

RITTER® EXAMINATION TABLES DESIGNED WITH YOU + YOUR PATIENTS IN MIND

What makes Ritter by Midmark the most purchased brand of examination tables in the industry? It's our willingness to research, study and listen so we can meet

your needs—and your budget. The Ritter 204 is the culmination of years of research in design, comfort and function in a manual examination table.

REMOVABLE TOP SYSTEM

The one-piece, upholstered top is designed to be easily removed without tools for a more thorough cleaning or a quick upholstery change.



The seamless design of the Ritter 204 upholstered top is as soft and comfortable

as it looks. And the seamless upholstery design reduces areas where dirt and bacteria can accumulate so clean-up can be easier, faster and more effective.

Stylish. Comfortable.

Easy to clean.

ONE-PIECE SEAMLESS DRAWERS

The Ritter 204 drawers are designed with rounded, seamless edges to help contain spilled fluids and improve infection control.



ENLARGED FOOTSTEP

We increased the size of the footstep to over 40% above the industry standard to help patients feel more secure when accessing and exiting the table.



SPACIOUS SIDE DRAWERS

Three side drawers provide point of care storage for gowns, drapes and other large items, as well as organized storage for smaller supplies and bulk items. Two front drawers provide easy access to items during lower body exams.



STORAGE ACCESS + VOLUME

Pass-through side drawers provide access to supplies from either side of the exam table and give you the flexibility to change your room layout or workflow.



PULLOUT STIRRUPS

Integrated in the exam table, the pullout stirrups provide leg support for patients during lower body examinations.





Upholstery Options



FlatProvides an additional working surface for exams and procedures where you need a little more space.



Plush, seamless upholstery designed to be more comfortable for the patient and easier to clean.



Soft TouchDesigned to be comfortable for the patient yet durable enough to stand up to the rigors of the clinical environment.

Accessories



Assist Handles
Designed to offer patients a
stable platform for accessing
the table and a sense of security
during ingress, egress and
positioning while on the table.



Ritter® 253 LED Exam Light
The 253 is designed with easy
access controls, high intensity
light and an adjustable focal
spot. Available in table or wall
mount and mobile options.



Adjustable Arm System
Designed to offer patients a stable platform in accessing the table and a sense of security during positioning. Arms rotate and lock into one of three positions and are easily removable.



Stainless Steel Treatment Pan Desigend for spaces where a stainless steel treatment pan is preferred.



Articulating Knee Crutches
The unique ball-joint design
allows for maximum range of
adjustment and can help ensure
optimum knee and leg support.



Standard Knee Crutches
Swivel to provide leg
support for extended
procedures where fatigue is
a factor or when the patient
lacks sufficient leg strength
or control.

Ritter® 203 Treatment Table

FEATURES

The economical Ritter 203 is constructed of steel from its frame to its legs. This metal-on-metal design helps give the 203 greater stability and structural integrity. Its flat-top design provides a stable, solid surface for patient manipulations or adjustments during examinations while extra padding and a built-in headrest pillow help keep patients comfortable. And the Mist colored base complements any of the Ritter upholstery colors.



OPTION

Drawers

Rounded drawers with color inserts are available on the right or left side of the table.



6 Manual Examination Tables 7

Ritter® 95 Treatment Table

FEATURES

The Ritter 95 is an economical solution that's built to function in a wide range of applications. It's ideal for new medical practices, industrial first aid areas, school nurses' stations or as an additional table wherever needed. The hardwood construction and H-brace frame design improve strength and stability while the 27" width and 72" length help accommodate a wider range of patients. (Paper roll holder included.)

OPTION

Vinyl Upholstered Shelf

Provides space at the point of care to place supplies or use as an additional worksurface.



8 Manual Examination Tables

SPECIFICATIONS

RITTER 204

Height with Upholstery: Soft Touch Top -32.13" (81.6 cm) Seamless Top -33.13" (84.2 cm)

Patient Weight Capacity: 500 lb (226.8 kg)

Electrical Receptacle: (115 VAC 5 A Max.) (230 VAC 2.5 A Max.) Duplex, hospital grade electrical receptacle on left side

Soft Touch Upholstery Top: 30.5" W x 56" L (77.5 cm x 142.5 cm)

Paper Roll Holder holds 21" x 4.75" Paper Roll (53.3 cm x 12 cm)

Length with Footrest extended: 73.75" (187.3 cm)

Seamless Upholstery Top: 28.5" W x 59.5" L (72.4 cm x 151.1 cm)

Paper Roll Holder holds 21" x 3.5" Paper Roll (53.3 cm x 8.9 cm)

Length with Footrest Extended: 76" (193 cm)

UltraFree Upholstery Top: 28.5" W x 59.5" L (72.4 cm x 151.1 cm)

OPTIONS

Soft Touch, Seamless, Seamless Ultrafree, Flat Stitched Upholstery without pillow, Flat Ultrafree Upholstery top without pillow

Pelvic tilt (1.5° or 7°)

Drawer heater (45 W)*

ACCESSORIES

Knee crutches

Stainless steel treatment pan

Assist handles

Adjustable arm system

Pass-through drawer stop (set of 3)

Thick footrest pad kit

ELECTRICAL

115 VAC, 50/60 Hz, 5.0 A (with receptacle)

115 VAC, 50/60 Hz, 5.5 A (with drawer heater)

230 VAC, 50/60 Hz, 2.5 A (with receptacle)

*Not available for sale in Europe.

RITTER 203

Patient Weight Capacity: 400 lb (182 kg)

Dimensions (top): 27" W x 74" L (68.6 cm x 188 cm)

Height: 31" (78.7 cm)

Paper Roll Holder: 21" (53.3 cm)

RITTER 95

Patient Weight Capacity: 350 lb (159 kg)

Dimensions (top): 27" W x 72" L (68.6 cm x 182.9 cm)

Height: 30" (76.2 cm)

Paper Roll Holder: 21" (53.3 cm)

ACCESSORIES

Vinyl Upholstered Shelf

COLOR OPTIONS

STANDARD UPHOLSTERY COLORS

Mist 813

Stone 814


Iron Ore 851

Oasis 843

Restful Path 858




Obsidian 857

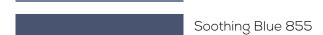







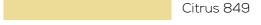











Latte 860

Deep Earth 854

ULTRAFREE UPHOLSTERY COLORS



11

10 Manual Examination Tables



CARB 93129.2 Phase 2 Compliant and TSCA Title VI Compliant

The color examples shown are the best representation of the original material. Actual color may vary slightly, We strongly recommend that you contact Midnark Customer Experience at 1.800.MIDMARK to request a sample before placing your order.

Authorized CE Representative Countries in the EEC should direct all questions, incidents and complaints to the Midmark Authorized CE representative listed below:

CEpartner4U, Esdoornlann 13, 3951 DB Maarn, The Netherlands Phone: +31.343.442.524 Website: cepartner4u.com

Midmark is an ISO 13485 and ISO 9001 Certified Company. Certain products are not included. See the complete list at: midmark.com/ISO

For more information, contact your Midmark dealer or call: 1.800.MIDMARK

Outside the USA call: 1.937.526.3662 or visit our website: midmark.com

© 2021 Midmark Corporation, Miamisburg, Ohio USA Products subject to improvement changes without notice. Litho in USA. **Ritter**®