


Anesthesia 101

Rebreathing Circuit

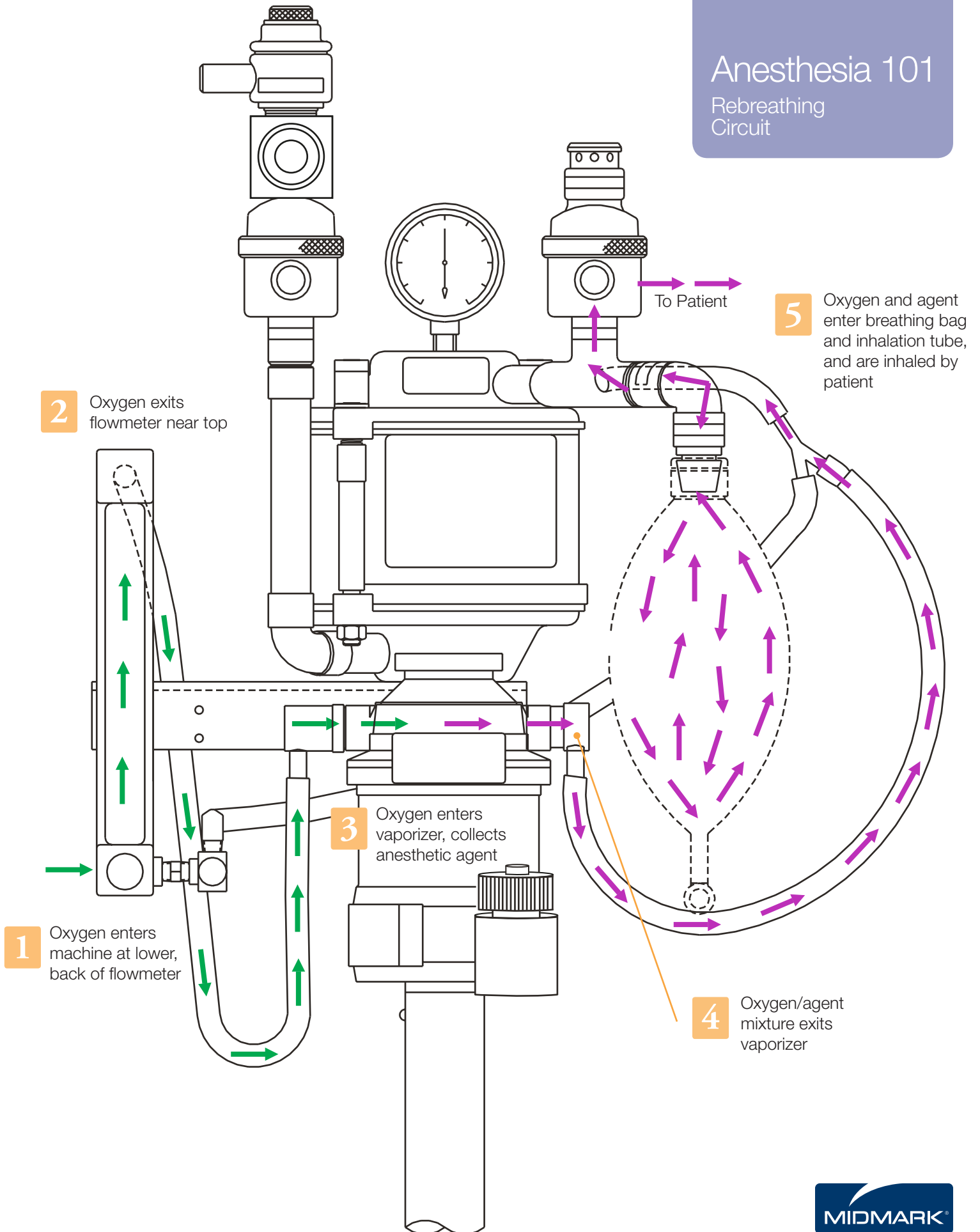

Anesthesia 101

Rebreathing Circuit

9

Small amount of waste gas exits APL valve

From Patient

6


Patient exhales - waste gas enters machine through exhalation valve

8

Gas with CO₂ removed enters breathing bag, inhalation tube and patient (recycled)

7

Waste gas enters CO₂ absorber canister, CO₂ removed

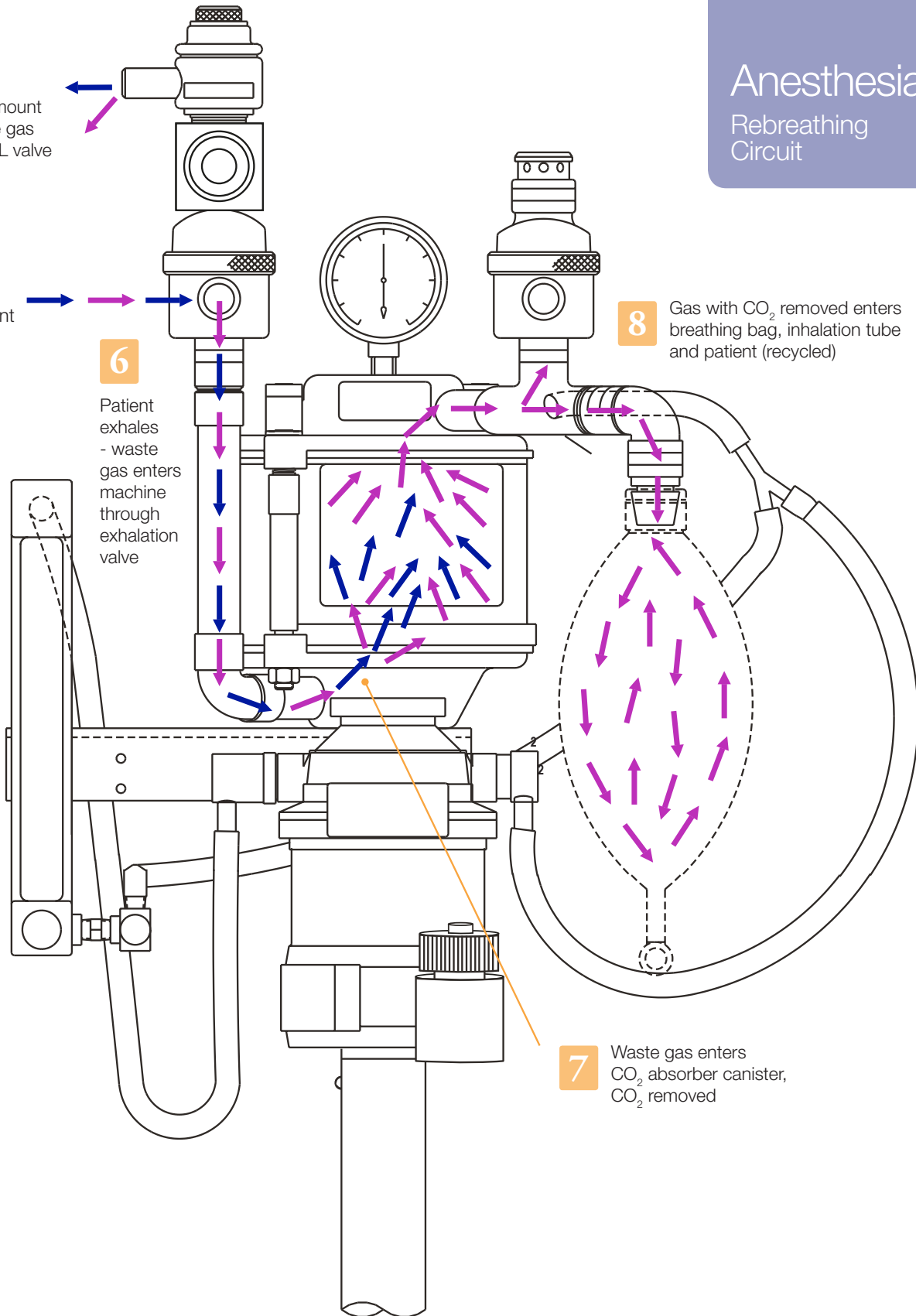

Better Patient Care. Better Business.™

   midmarkanimalhealth.com

Manufactured and/or distributed by Midmark Corporation, Versailles, OH.

© 2005 Midmark Corporation • Products subject to improvement changes without notice • 007-1073-00 Rev. B1 (1/16)
Midmark is an ISO 13485 Certified Company.



Because we care.