

Synthesis® Fixed-Height Wet Treatment Tables

Designed
to **Support**
the Way
You Work




Synthesis Fixed-Height Wet Tables are built for the realities of veterinary care—whether you're treating, cleaning or prepping patients. Blending thoughtful design with lasting durability, they help your team stay focused, efficient and comfortable.

- Choose a peninsula-style wall or chase mount, or freestanding island configuration to fit your space and workflow.
- Accommodate different treatment needs with 5" deep or 3" shallow tub options.
- Enhance patient comfort with flat, black, powder-coated grates that provide greater surface area for sensitive paws.
- Personalize your space with more than 20 cabinet color options designed to complement your facility's style.
- Maintain a quiet, low-stress environment with soft-close hinges and slides.
- Streamline workflow with a multi-grate design—use only the grates you need to improve flexibility and reduce cleaning time.
- Simplify cleaning with a stainless steel tub and integrated faucet and sprayer.
- Support infection control protocols with optional antimicrobial pulls and surfaces.



NOTE: Catalog-ordered Synthesis Wet Tables ship standard with black flat grates. When ordering through your sales representative, you may choose either black flat or stainless steel wire grates.

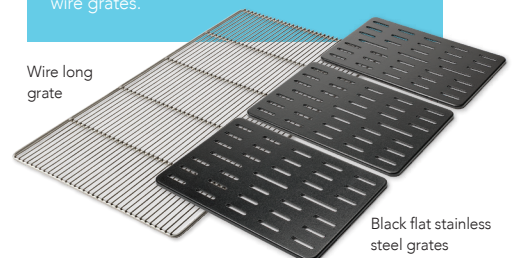


Wire long grate


Black flat stainless steel grates

SYNTHESIS® FIXED-HEIGHT WET TREATMENT TABLE SPECIFICATIONS

STANDARD Wet Table and Plumbing Cutout Dimensions

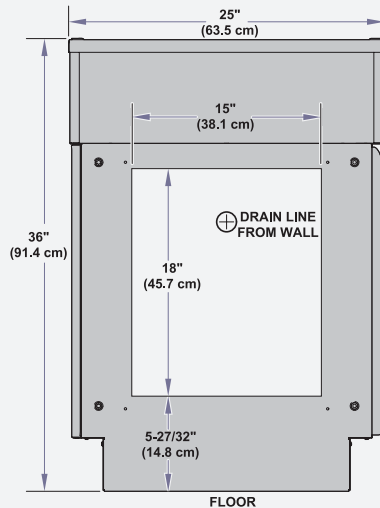

SHALLOW Wet Table and Plumbing Cutout Dimensions

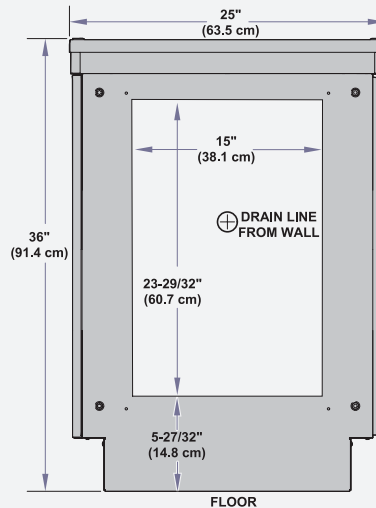

TECHNICAL SPECIFICATIONS

Height	36"
Length	48" or 60" (includes 18" overhang)
Width	25"
Weight	250 lb (48" table) 275 lb (60" table)
Weight Capacity	300 lb
Frame	Thermofoil over electrostatic, powder-coated painted steel (woodgrain and metallic colors); electrostatic, powder-coated painted steel (solid colors)
Doors and Drawers	Sealed ¾" medium-density fiberboard with PVC thermofoil and electrostatic, powder-coated paint
Configurations	– Tub: 5" deep or shallow 3" – Left or right-hand layout – Multiple cabinet door and drawer options – Available with plumbing cutout or solid back
Grate Style Options	– Black Flat Stainless Steel Grates (standard with catalog orders; 2 grates with 48" table, 3 grates with 60" table) – Stainless Steel Wire Long Grate
Accessories	Faucet and sprayer included



- 01** Durable 18-gauge steel frame protects against rust, while sealed cabinet panels resist moisture and chemical damage.
- 02** Seamless polystyrene drawers with rounded corners and one-piece construction allow drawers to contain spills.
- 03** Modular drawer and door layouts make it easier to organize supplies.

ORDERING GUIDE

For quick ordering, explore the many pre-designed table options in the Midmark Equipment Catalog, or contact your Midmark sales representative to custom-configure the right solution for your space and care.



Scan QR code to view the Midmark Equipment Catalog.