


Dental/Surgery Flat-Top Table

Designed
for patient +
care team
comfort + safety



The Midmark Dental/Surgery Flat-Top table offers unencumbered access to the patient during procedures such as dentistry. Enhanced ergonomic and hygienic features promote an optimal care environment for patients and caregivers.

- Durable metal grate and built-in drainage basin help ensure fluids efficiently drain away from patient.
- Raised marine edge prevents fluids from dripping onto the floor.
- Removable 3-quart fluid receptacle is easy to clean and disinfect.
- 50" and 60" table lengths accommodate a variety of procedure spaces.
- Adjustable table height and tilt support proper ergonomics in both electric and hydraulic models.
- Optional thermal heating helps maintain the patient's body temperature.
- Optional heavy-duty casters allow you to effortlessly move the table to where it is needed.

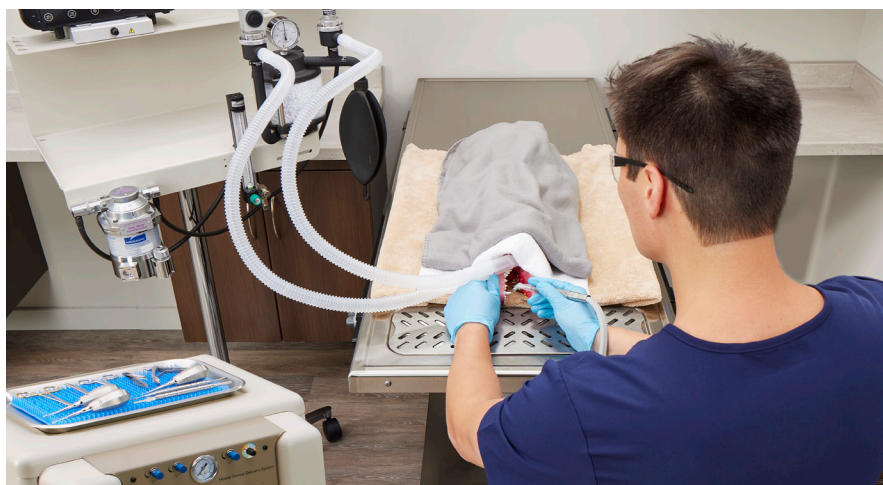

CREATE YOUR OWN BEST-PRACTICE DENTAL SUITE

Combine the Midmark Dental/Surgery Flat-Top Table with our dental delivery system, digital dental radiography, anesthesia, multiparameter monitoring and Synthesis® Cabinetry.

MIDMARK DENTAL/SURGERY FLAT-TOP TABLE SPECIFICATIONS



- Length:** 50" and 60" models
- Width:** 22" on both table models
- Height:** Adjustable from 29 ½" – 41 ½" on electric model and 29 ½" – 35" on hydraulic model
- Fluid Receptacle:** 96 oz at 0 tilt
- Optional:** Heated or Non-Heated Top

ORDERING GUIDE

Item Number	Description:
100-2091-00	Dental/Surgery Flat-Top Table, 60", Electric, Heated
100-2092-00	Dental/Surgery Flat-Top Table, 60", Electric, No Heat
100-2093-00	Dental/Surgery Flat-Top Table, 50", Electric, Heated
100-2094-00	Dental/Surgery Flat-Top Table, 50", Electric, No Heat
100-2095-00	Dental/Surgery Flat-Top Table, 60", Hydraulic, Heated
100-2096-00	Dental/Surgery Flat-Top Table, 60", Hydraulic, No Heat
100-2097-00	Dental/Surgery Flat-Top Table, 50", Hydraulic, Heated
100-2098-00	Dental/Surgery Flat-Top Table, 50", Hydraulic, No Heat

