

# How to Update the Midmark PDF Converter driver

## Problem

---


This knowledge article details the process for updating the Driver of the Midmark PDF Converter from Amyuni Document Converter 400 to Amyuni Document Converter 450

## Solution

---

### Step 1


- Open the .cab file by double clicking it
- Highlight all the contents of the .cab file (see screenshot below) and copy them

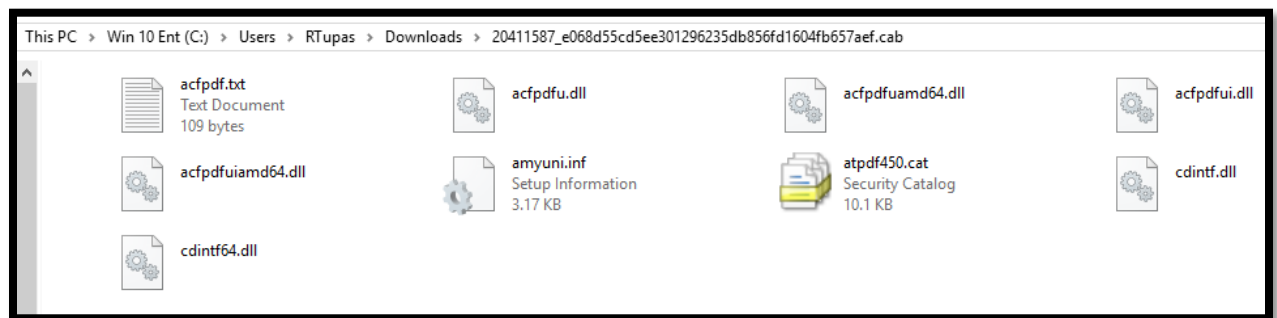


- Paste them to a specific folder of your choice which you have full access and rights

### Step 2


- Launch Control Panel
- Select Devices and Printers from the item list
- Locate the Midmark PDF Converter in the Printers region
- Right click on the Midmark PDF Converter and select Printer Properties

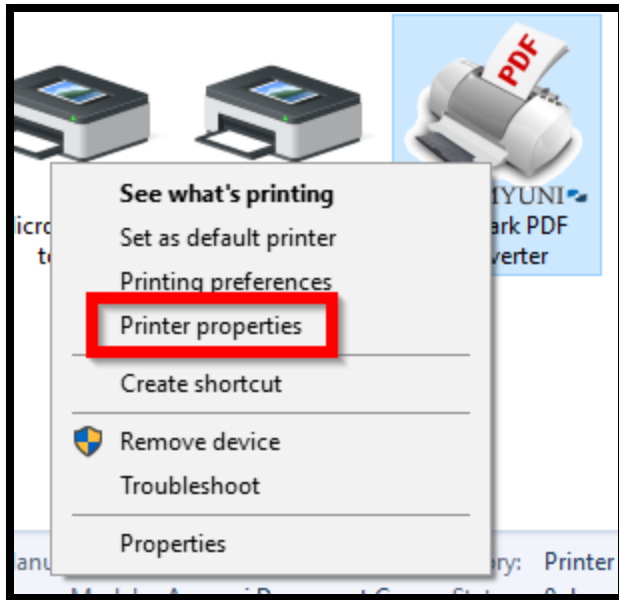


- Select the Advanced tab

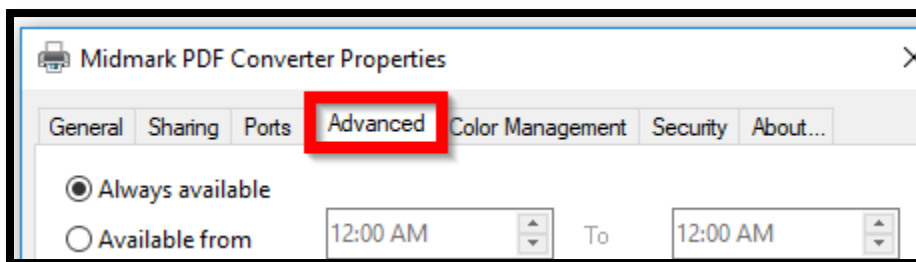


- In the driver region click the New Driver button
- The Add Printer Driver Wizard will open then just click on the Next button
- The next window will ask if you have the driver under the Windows Update or if you Have Disk
- Select/click the Have Disk button
- Browse to where you extracted the contents of the downloaded .cab file
- Select the Amyuni Document Converter 450 from the list

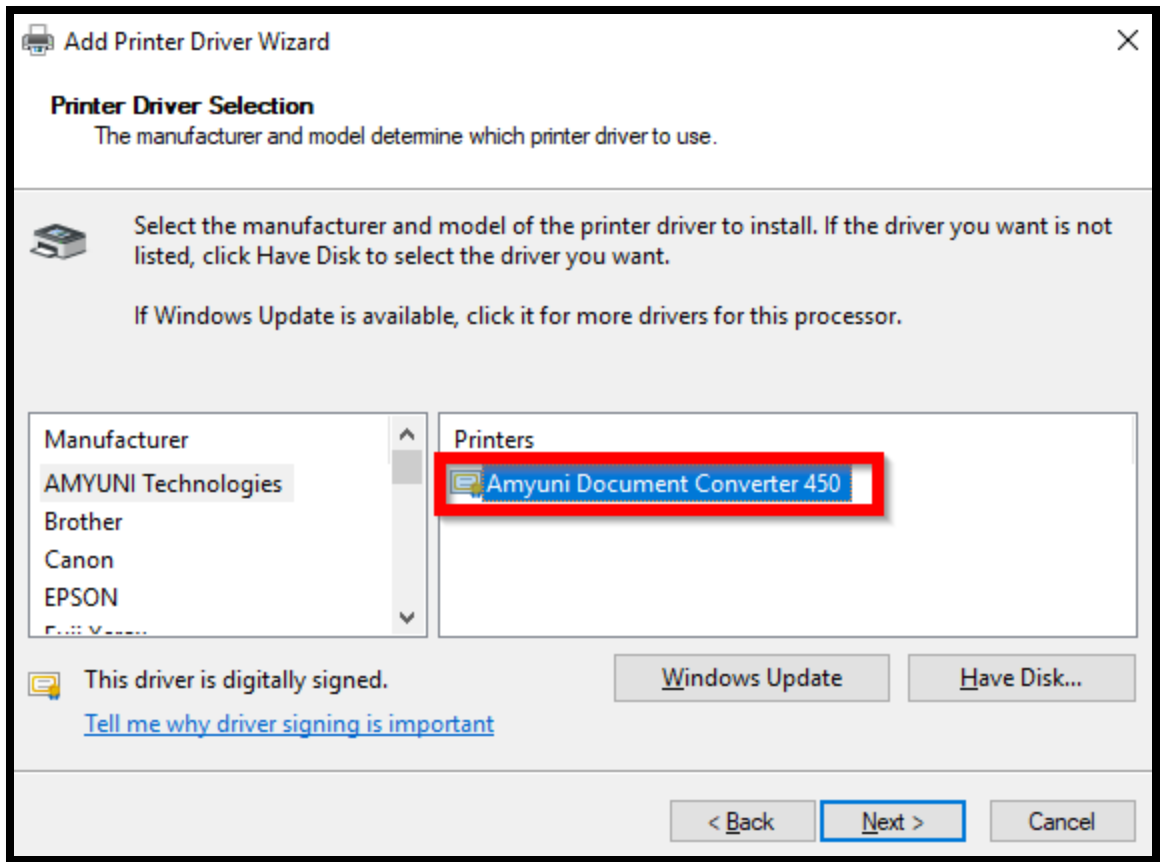

- Select/click the Next button
- Select/click the Finish button
- You should now see the selected driver as Amyuni Document Converter 450

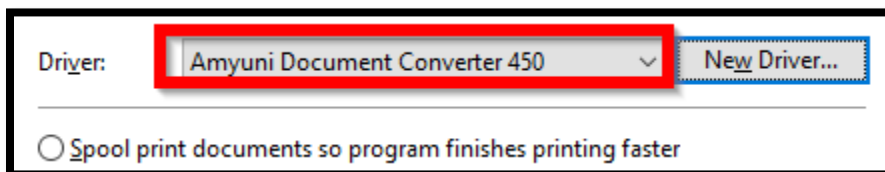

### Step 3

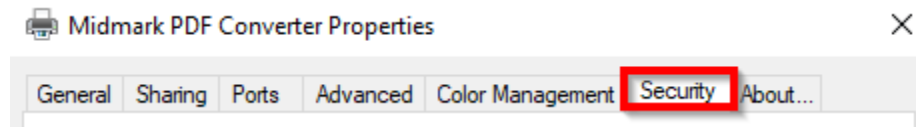

- Select/Click the Security tab
- Under Group or user names, make sure all listed groups/user names have full permissions (Allow)
- Once all group/user names have full permissions, select/click the OK button

