

# Ritter®

MEDICAL

## Direct-To-Drain System




# Say goodbye to your condensate tanks.

The steam and condensation produced by your sterilizer can be time-consuming and messy to empty from condensate tanks. It can even melt plumbing if not cooled sufficiently before pouring down the drain.

There is a better solution. The VistaCool™ Direct-to-Drain System for Autoclave Wastewater is designed to reduce the temperature of sterilizer wastewater—eliminating condensate bottles and tanks!




## SINGLE UNIT

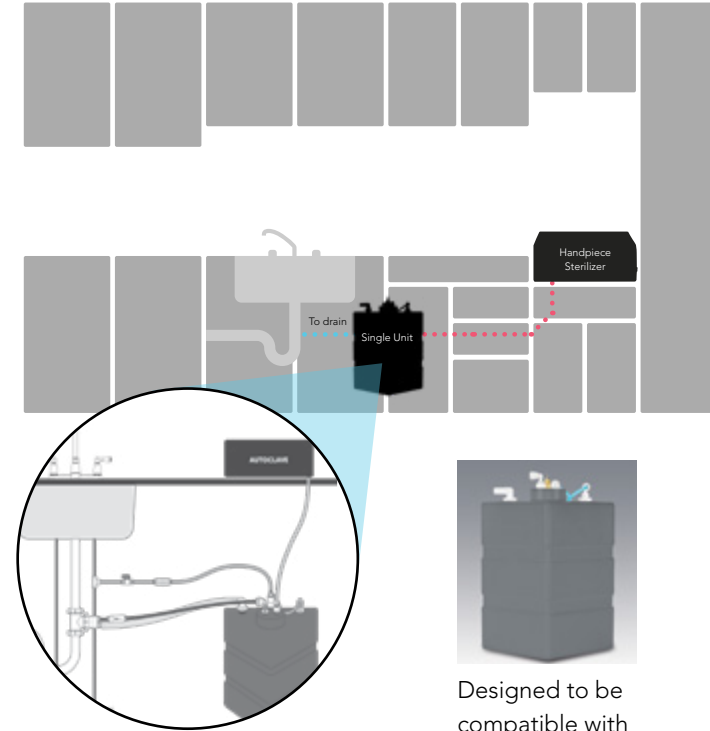

Designed to be compatible with the Ritter M3 Part No. 9A586001

## DOUBLE UNIT

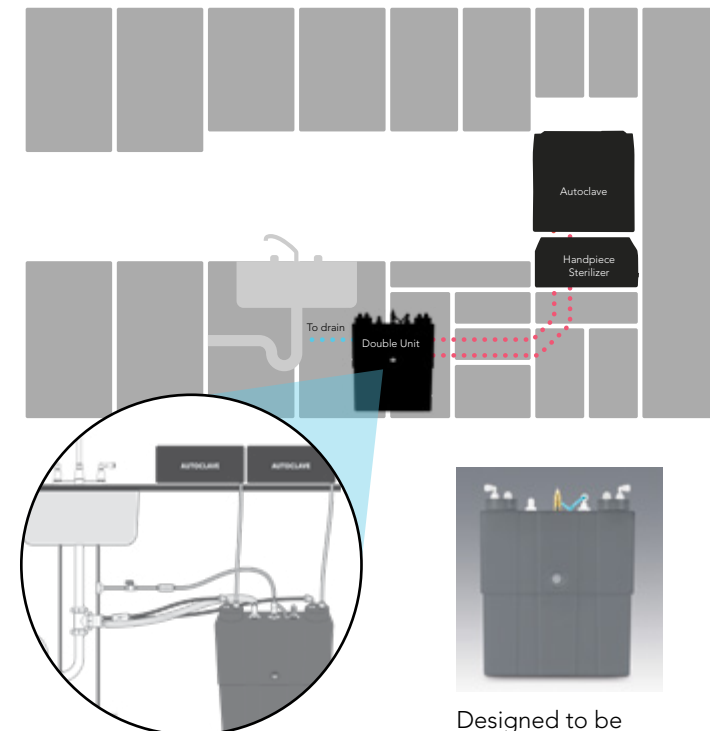

Designed to be compatible with any two current Ritter sterilizers Part No. 9A586002

*VistaCool is the only direct-to-drain product designed and tested to work with Ritter automatic sterilizers.*

# VISTACOOOL™ IS MORE THAN A PRODUCT, IT'S A SOLUTION.

Eliminating the need to manually empty condensate tanks, VistaCool will allow you and your staff to spend more time with patients and less time tending to office equipment. This safe and reliable solution is designed to enhance office efficiency. Less work for staff means more time to care for patients.

## PROTECTS STAFF

- Eliminates hazardous handling of hot condensate tanks
- Non-electric, self-monitoring and self-regulating system saves staff time with limited user maintenance required.

## PROTECTS DRAINS

- Designed to prevent sterilizer wastewater from melting plumbing
- Built-in backflow prevention meets or exceeds local plumbing codes

## PROTECTS CABINETS

- Helps prevent mold and mildew caused by spills and exhaust from condensate tanks
- Helps protect cabinets from rot and rust by sending cooled wastewater directly and safely to the drain

Designed to be compatible with:  
Ritter M3 Steam Sterilizer  
Ritter M9 and M11 Steam Sterilizers\*

\*Required adapter kit included. Compatible with the double unit only.

# SPECIFICATIONS

	Single	Double
Compatible Sterilizer	M3	M3, M9 and M11
Tank material:	Polyethylene - Gray	Polyethylene - Gray
Tank dimensions w/ fittings:	7.6" W x 7.6" L x 13.75" H (19.3 cm x 19.3 cm x 34.9 cm)	4.0" W x 13.0" L x 17.0" H (10.2 cm x 33.1 cm x 43.2 cm)
Tank volume:	1.50 gal (5.70 liters)	3.00 gal (11.40 liters)
Tank weight w/ water + fittings:	15 lb (6.8 kg)	28 lb (12.7 kg)
Approximate shipping weight:	7 lb (3.2 kg)	10 lb (4.5 kg)
Coolant water connection:	Supply fittings for 1/4" O.D. feed tubing included. To be provided by installer.	Supply fittings for 1/4" O.D. feed tubing included. To be provided by installer.
Max. coolant water pressure:	100 psi (6.9 bar)	100 psi (6.9 bar)
Max. altitude (If using air gap):	6,561 ft (2,000 m)	6,561 ft (2,000 m)

Accessories		
9A586004 Air Gap Kit	(needed only if required by local code)	(needed only if required by local code)
9A586005	Single VistaCool adjustable tank stand	N/A
9A586006	N/A	Double VistaCool adjustable tank stand
*9A586003	N/A	Optional M9/M11 Adapter Kit

\*Required adapter kit (9A586003) not included. M9 and M11 Steam Sterilizers are compatible with the double unit only.

## VistaCool™

Midmark is an ISO 13485 and ISO 9001 Certified Company.<sup>†</sup>  
 For more information, contact your Midmark dealer  
 or call: 1.800.MIDMARK  
 Outside the USA call: 1.937.526.3662 or visit our  
 website at [midmark.com](http://midmark.com)

© 2019 Midmark Corporation, Miamisburg, Ohio USA  
 Manufactured by Crosstex International, Inc.,  
 for distribution by Midmark Corporation, Versailles, OH  
 VistaCool™ Direct-to-Drain System for Autoclave Wastewater is a trademark  
 of Crosstex International, Inc., a Cantel Medical Company, Hauppauge, NY  
 Products subject to improvement changes without notice.  
 Litho in USA  
 007-10159-00 Rev. A1 (7/19)

<sup>†</sup>Certain products not included.  
 See the complete list at: [midmark.com/ISO](http://midmark.com/ISO)

# Ritter®