


Monitor Mount

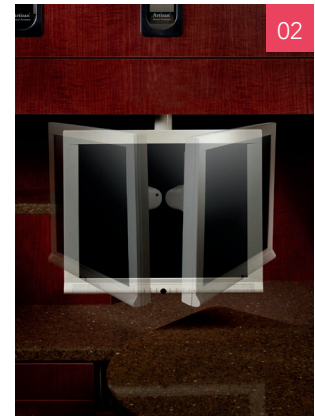

Improved efficiency. Better ergonomics. More workspace.

The Midmark sliding track monitor mount is designed to assist dental clinicians with positioning a monitor on a Synthesis® and Artizan® Expressions rear cabinet, providing an efficient and ergonomic method of viewing media such as office schedules, patient charts and digital x-rays. The mount allows the monitor to be suspended above the countertop, providing more space to work and extra room for essential supplies and ancillary devices.

Specifications

- Maximum monitor weight: 15 lb
- Accommodates up to 24" monitor
- Vertical adjustment: 5"
- Swivel range: 180°
- Tilt range: 15° forward, -15° backward
- Track length: 22"

Features:

- 01 Side-to-side and wrap-around sliding track allows mount to slide from doctor's side to assistant's side and also wrap around the cabinet to be viewed outside the operator.
- 02 Left-right pivoting of monitor assists dental clinicians in viewing media from a variety of angles.
- 03 Tilt feature helps reduce glare when viewing from a standing or seated position.
- 04 Height Adjustment mechanism allows the monitor to be adjusted to avoid support trays and taller devices.

Additional Features:

- Easy positioning of monitor
- Accommodates both 100/75 VESA mounting patterns
- Finished back cover
- Cable management system
- Cable options available
- 2-year limited warranty